


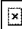
From: Separation Science <elearning.solutions@sepscience.com>
Sent: Friday, June 29, 2012 1:09 PM
To: Hanchett, James (DPH)
Subject: Today in Separation Science - Latest Issue of HPLC Solutions



Web Version Forward Unsubscribe

John Dolan's
Update
HPLC Solutions

Separation Science e-Learning


Agilent Technologies is pleased to announce the release of the latest issue of the HPLC Solutions newsletter.

Issue 97
Late Elution
A reader asks, "I have a nice isocratic method for compound A, which elutes at about 8 min, so I can inject a new sample every 10 min. Occasionally, however, my samples contain another compound B that comes out at about 40 min. I have no way of knowing in advance if B will be present or not. Some people in the lab are telling me that I need to wait 40 min for every sample to be sure that this peak is not present before I can inject again. This seems like a lot of wasted time. Can't I just reinject the sample when B shows up in the wrong place in its chromatogram? After all, the B peak is pretty broad and easy to recognize."
[Click for PDF>>](#) [View Arch!](#)
[Ive>>](#)

Issue 96
Premix
A reader asked me when it is appropriate or necessary to premix mobile phases when on-line mixing is available. For example, if a gradient is run from 95/5 buffer/organic (5% B) to 5/95 buffer/organic (95% B), should one bottle contain 95% buffer and the other 5% buffer in the organic solvent, or should one bottle contain buffer and the other the organic? In either case, the HPLC system can be programmed to give exactly the same mobile phase composition.
[Click for PDF>>](#) [View Archive>>](#)

Issue 95
UV Conditions
I recently had an "Ask the Doctor" question that went something like this: "I need to report impurities for my product at levels of 0.1-0.5% of the main ingredient. Occasionally when I show my folks my chromatograms, I'm accused of making the sample too concentrated, because the product peak is >2000 mAU. They tell me that the detector is saturated, so the area-% for the impurities will not be accurate. Can you comment on this?"
[Click for PDF>>](#) [View Archive>>](#)

ASK THE DOCTOR
If you have an analytical question for John Dolan...
[CLICK HERE >>](#)

FEATURED APPLICATIONS

Single-run assay and impurity testing of fixed-dose combination drugs using the Agilent 1200 Infinity Series High Dynamic Range Diode Array Detector Solution
Company: Agilent Technologies
[Click to read>>](#)

Analysis of 'Broad Spectrum' UVA and UVB Components in Sun Care Products for Compliance with New FDA Regulations
Company: Perkin Elmer
[Click to read>>](#)

Extraction of Malachite Green, Crystal Violet and their Leuco Metabolites From Salmon Using QuEChERS and EVOLUTE CX for LC-MS-MS Analysis
Company: Biotage
[Click to read>>](#)

Analysis of Intact Proteins on a Thermo Scientific Accucore 150-C4 150 Å Pore Diameter HPLC Column
Company: Thermo Scientific
[Click to read>>](#)

Analysis of Pioglitazone Hydrochloride
Company: GL Sciences
[Click to read>>](#)

Published by Eclipse Business Media Ltd
Frederick House | Princes Court | Beam Heath Way | Nantwich | Cheshire CW5 6PQ | United Kingdom
20 Maxwell Road | #06-17 Maxwell House | Singapore 069113

Copyright © 2012 Eclipse Business Media Ltd. All rights reserved.

To sponsor or advertise in this e-newsletter contact dean.grimes@sepscience.com
For editorial contributions, contact david.hills@sepscience.com

This message was sent from Separation Science to james.hanchett@state.ma.us. It was sent from: Eclipse Business Media Ltd, Frederick House, Princes Court, Beam Heath Way, Nantwich, Cheshire CW5 6PQ, United Kingdom. You can modify/update your subscription via the link below.



[Unsubscribe](#) | [Forward To a Friend](#)